
Withings Thermo

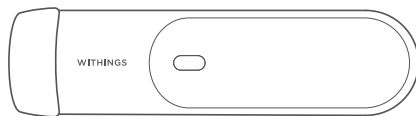QUICK INSTALLATION GUIDE

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Read this manual before use

EN | SW | FI

THANK YOU FOR CHOOSING WITHINGS THERMO

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Behöver du hjälp? | Tarvitsetko apua?

➔ withings.com/support/thermo

This guide is available in additional languages at: www.withings.com/guides

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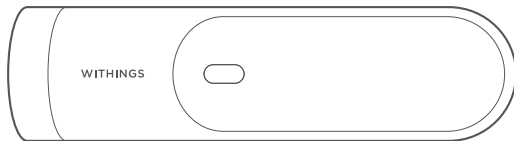
How to install your Withings Thermo

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4-5

1

Visit thermo.withings.com from your smartphone or tablet, and follow the instructions to install the Withings Thermo app.

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
2

In the Thermo app, follow the instructions to install Thermo.

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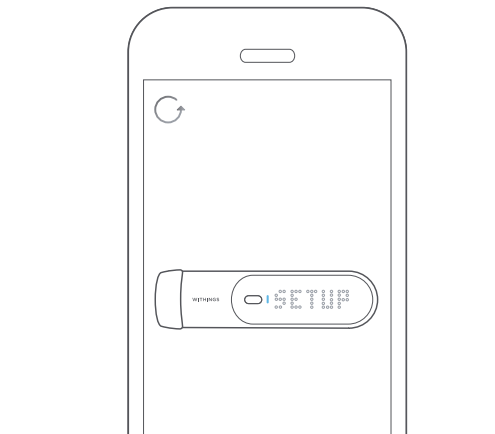1

Remove the green protective cap.

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
i Clear any sweat or hair from the forehead before taking a measurement.

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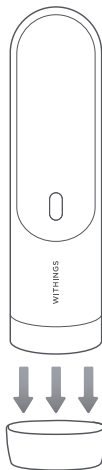2

Press the button to turn the thermometer on.

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
Clear any sweat or hair from the forehead before taking a measurement.

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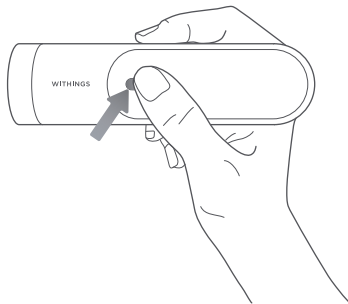
3

Starting from the middle of the forehead, slowly scan straight across to the hairline, up to 1cm over the skin.

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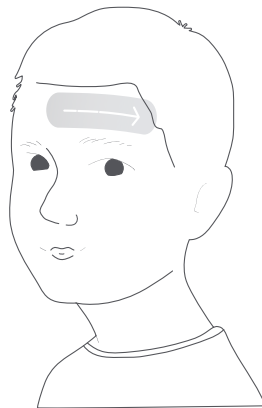
4

Thermo will vibrate twice when the measurement is complete. The value of the temperature is displayed on the device. The unit of this value is either °F or °C. You can change this unit in the Thermo app settings. If you change the unit, make sure to update your device in the Thermo app settings.

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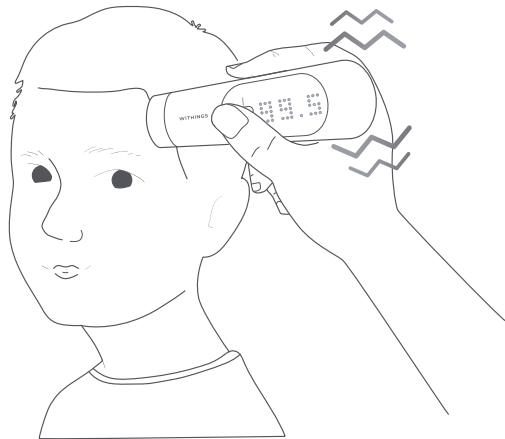
5

To assign the temperature to a user, scroll through names on the touch-sensitive zone (1). Press the button to assign the reading (2).

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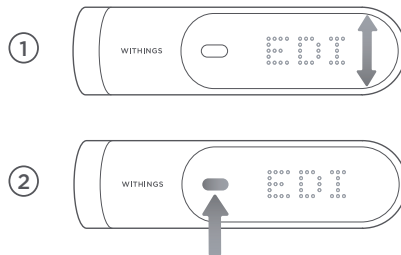
6

A colored LED light indicates the fever level according to the user's age.

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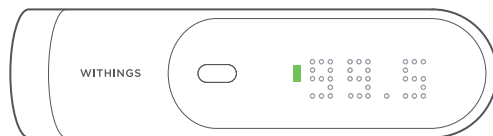T (°C)

Important advice

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The person being measured must be indoors for 15 minutes prior to taking a reading. | The person being measured must have dry skin. If needed, dry skin with a towel prior to taking a reading. | The person being measured and the thermometer should be in the same ambient temperature for at least 15 minutes prior to taking a reading. | The person being measured should not have exercised within 15 minutes prior to taking a reading. | Don't take a measurement while or immediately after nursing a baby.

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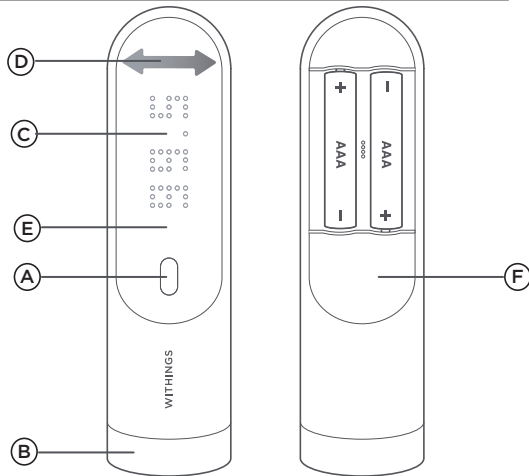
Warranty

Garanti | Takuu

UK Withings Two (2) Years commercial guarantee - Withings Thermo Withings SA, 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, France («Withings») warrants the Withings-branded hardware product («Withings») against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of TWO (2) YEARS from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

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EN Withings SA hereby declares that the device Withings Thermo is in conformity with the essential requirements and other relevant requirements of Directive 2014/53/EU. The declaration of conformity can be found at: withings.com/compliance

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The crossed-out wheeled bin symbol on your product, battery and/or accessories means that the item is classified as electrical and electronic equipment.

QIG_SCT01_64

WITHINGS

Withings Thermo | Smart temporal thermometer - Model: SCT01 - V1.1

© 2020 Withings SA. All rights reserved.